

CONTINUE TO HEAL

THE LATEST MEDICAL TREADMILL FROM WOODWAY

Experience a lower RPE (Rated Perceived Exertion), less impact and less muscle fatigue with the CONTINUUM. With a low profile 8" step height and an absolute true zero start speed, patients will appreciate the safety and usability that the most comfortable and accessible rehabilitation platform has to offer.




CONTINUUM



- Standard Feature
- Optional Feature

		CONTINUUM
USER INTERFACE	Coded/Non-Coded Polar Heart Rate Pick-up	●
	ANT+ Heart Rate Pick-up	●
	USB Power Supply	
	Touch Grip Heart Rate Pick-up	
	LED Display Board	●
	Personal Trainer Display Board	○
	RS 232 Serial Interface	●
	CSAFE Communication Compatible	
	Wired Remote	●
	MED Pro Software	○
RFID System	○	
SPECIFICATIONS	Running Surface Dimensions	22" W x 54" L (55 x 137 cm)
	User Weight Capacity	Run: 400 lb. / Walk: 800 lb. (4 mph max)
	Belt Type	47 Individual Slats
	Drive System	80 Precision Ball Bearings with 8 Roller Guides (4 mm lateral tolerance)
	Running Surface	Vulcanized Rubber (38-43 shore hardness)
	Drive Motor	2 hp Continuous (5 hp peak) Brushless Servo
	Unit Weight	460 lb. (Shipping Weight 560 lb.)
	Width	39" (99 cm)
	Length	71" (180 cm)
	Height	60" (152 cm)
	0-25% Incline	●
	0-10 MPH (0-16 km/h)	●
	Reverse	●
	115 Vac 20 Amp	●
	208/230 Vac 20 Amp	○
	Hospital Grade Circuitry	●
	Warranty	5 Year Drive, Motor, & Belt 4 Year All Components 1 Year Labor

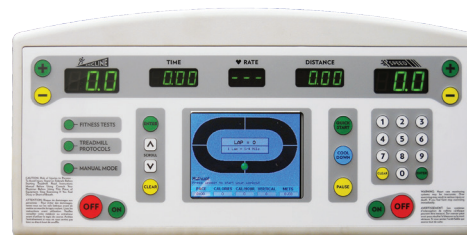

*Personal Trainer Display Option

CONVENIENT USER CONSOLE

- Comfortable handrail control of speed, elevation, and stop
- Multiple LED readouts monitoring speed, incline, distance, time, and heart rate
- CSAFE fitness communications compatible
- Hand-held remote (speed, incline & stop)
- Pre-loaded & customizable protocols & programs (Balke, Bruce, Gerkin, Cooper)



WOODWAY SLAT BELT RUNNING SURFACE

- Individual slats with vulcanized rubber surface
- Zero stretch & Zero slip toothed, wire-reinforced lateral belts
- Near frictionless precision ball bearing rail guide

SAFETY & RELIABILITY

- Emergency stop magnet and safety lanyard
- Hospital grade circuitry
- Polar monitor circuitry – integrated heart rate monitoring
- Advanced Slat Belt shock absorption
- Life expectancy of running surface 150,000 to 200,000 miles
- Extended warranty and factory renewal program

*Treadmill specifications subject to change without written notice.
Warranties outside the U.S. may vary - Please contact your dealer for details.*

*Full Compliance with the Legal Requirements for
Medical Products in the European Economic Area
(CE marking)*

WOODWAY Continuum treadmills comply with the following:

- Medical Devices Directive (directive 93/42/EEC)*
- Medical Standard EN 60601-1, 3rd edition
- Medical Standard UL 60601-1 (USA)
- Medical Standard UL CSA C22.2 No 601.1 (Canada)

**FOR ADDITIONAL INFORMATION ABOUT THIS PRODUCT
VISIT WWW.WOODWAY.COM OR CALL 800.WOODWAY TO TALK TO A REPRESENTATIVE.**